

# CS-Studio setup for areaDetector


Wayne Lewis

Osprey DCS

2018-03-08


# CS-Studio OPI Editor Perspective

- Open Perspective menu

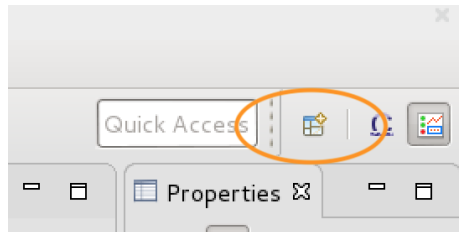

# CS-Studio OPI Editor Perspective

- Select OPI Editor
- OK

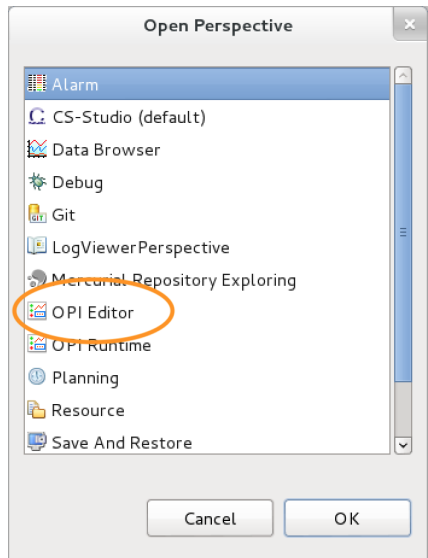

# Getting BOY displays

Get display repository. Clone location is arbitrary, but remember it.


[https://github.com/waynelewis/201803\\_training\\_opi](https://github.com/waynelewis/201803_training_opi)

- `cd ~/build-epics`
- `mkdir git`
- `cd git`
- `git clone`  
`https://github.com/waynelewis/201803_training_opi.git`

# Import project into CS-Studio


## Import into CS-Studio

- Launch CS-Studio:
  - ▶ `cd ~/cs-studio`
  - ▶ `./cs-studio-with-pref`
- Import new files:
  - ▶ File -> Import...

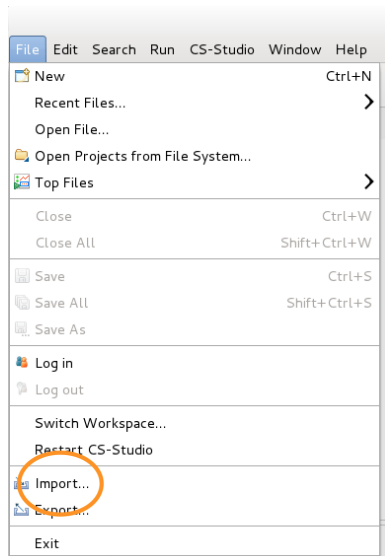


# Import project into CS-Studio

- General -> Existing Projects into Workspace -> Next>

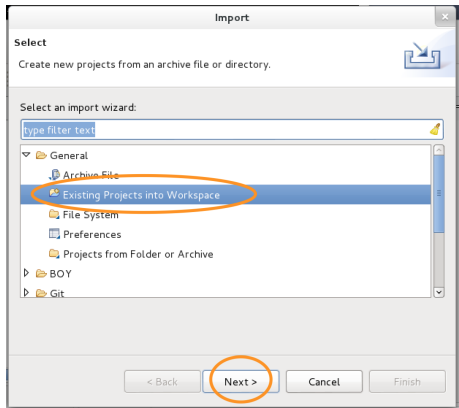


# Import project into CS-Studio

- Select root directory: -> Browse

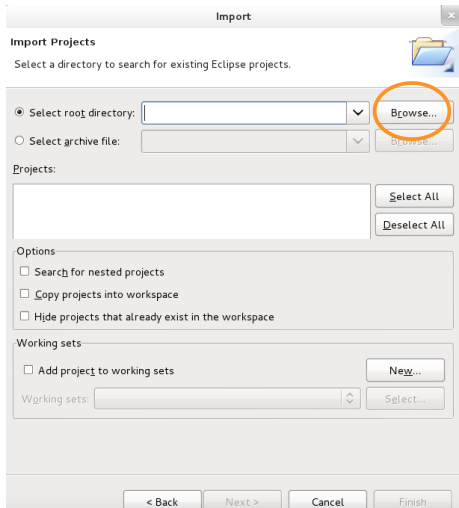

# Import project into CS-Studio


- Navigate to 201803\_training\_opi directory cloned earlier
- OK

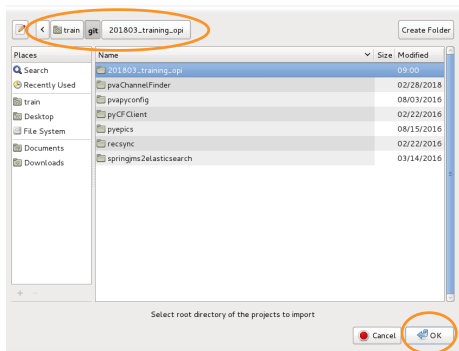




# Import project into CS-Studio

- Check correct directory location in 'Select root directory'
- Select 201803\_training under Projects: heading
- Finish

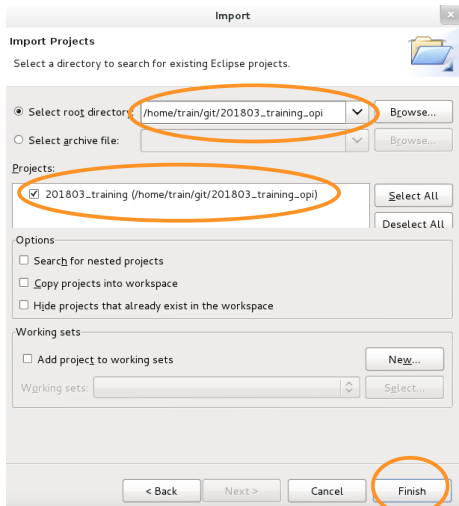


# Import project into CS-Studio

- Successful import

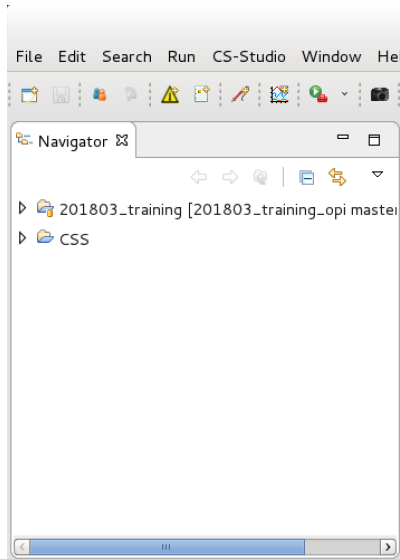


# Open OPI display

- Navigator -> 201803\_training
- Right-click -> Context menu
- Open With -> OPI Editor

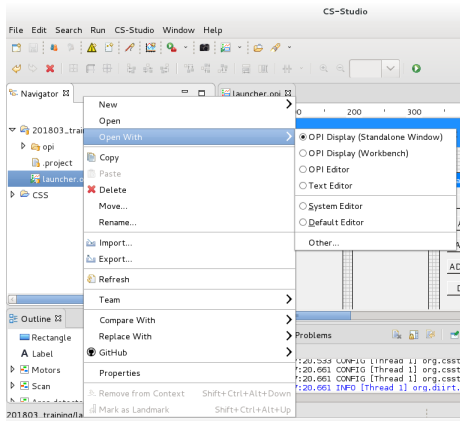

# Start OPI display

- Ctrl-G or Run OPI button

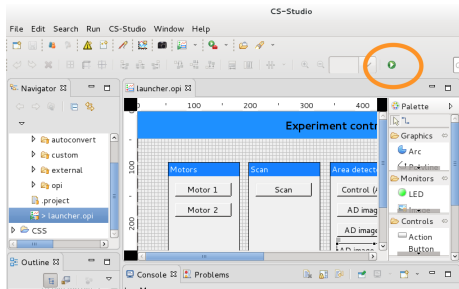

# Setting Max Array Bytes in CS-Studio

Default CS-Studio Max Array Bytes  
= 16384 - not enough for most  
images.

To set CA Max Array Bytes

- Edit -> Preferences
- CS-Studio Core -> Data Sources  
-> Channel Access

